

Direct, Rapid Analysis of Orange Peel using the ASAP Ion Source



Objective


To directly and rapidly analyze contaminants from orange peel with high mass accuracy and resolution using an ambient desorption ionization technique, Atmospheric Solids Analysis Probe (ASAP), and an Agilent TOF or Q-TOF instrument.

Background

ASAP uses heated nitrogen to vaporize the sample and a corona discharge for ionization. The ion source is a modified APCI source with a corona needle and a probe for sample introduction. This source only takes a few minutes to set up and run.

This technique is fast, easy to use and does not require sample preparation. Using a high resolution instrument such as a TOF or a Q-TOF, it is possible to acquire high resolution accurate mass MS and MS/MS data in less than a minute.

This note describes the analysis of fungicides such as imazalil and thiabendazole from orange peel in less than 30 seconds without sample preparation using ASAP and an Agilent Q-TOF instrument.

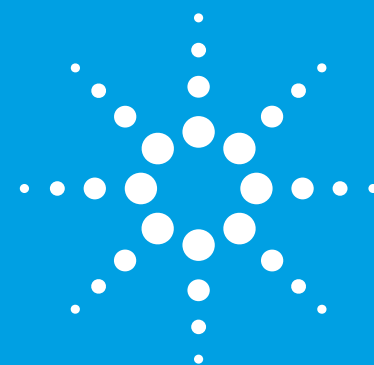

Compounds

- Imazalil
- Thiabendazole

Key Benefits

- No sample preparation is needed with ASAP analysis
- Complex samples can be analyzed in less than 30 seconds using ASAP
- ASAP is fully compatible with Agilent TOF and Q-TOF instruments
- Typical mass accuracies for MS and MS/MS data are within 3 ppm



The Approach

The spectra shown below were acquired in less than 30 seconds following exposure of the sample in the ionization source. For ASAP analysis, a very small amount of sample can be placed in a sampling glass capillary. The capillary is then mounted onto a solids probe and introduced into the APCI source.


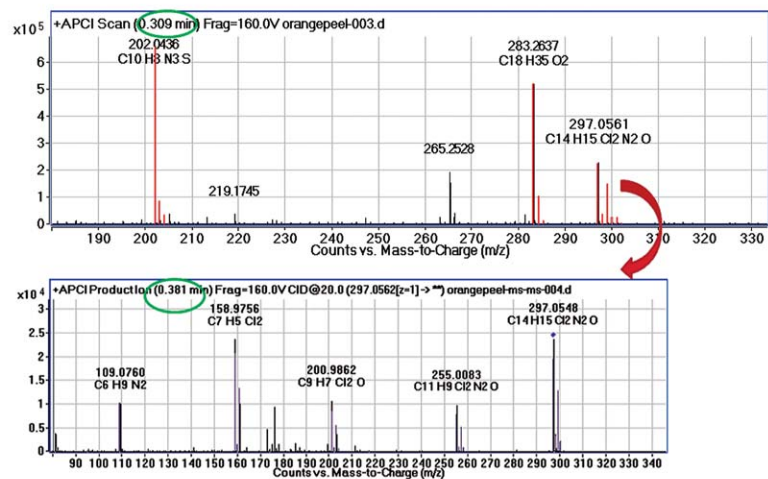


Figure 1. Orange peel scraped and analyzed with ASAP and an Agilent Q-TOF. The MS spectrum is shown in the top panel. Mass accuracy is within 2 ppm. The MS/MS spectrum of imazalil from orange peel is shown on the bottom panel. Mass accuracy is within 3 ppm. The data acquisition time of 0.3 min is circled in green in the top portion of the chromatogram.

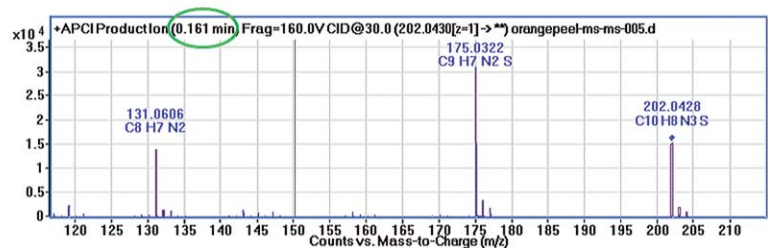

MS Formula Results: • Product Ion (0.161 min) (202.0430 > *)

m/z	Ion	Formula	Abundance
131.0606	M ⁺	C8H7N2	13963.4
175.0322	M ⁺	C9H7N2S	31257.5
202.0428	M ⁺	C10H8N3S	15440.1

m/z	Ion	Formula	Abundance
131.0606	C8H7N2	C8H7N2	13963.4
175.0322	C9H7N2S	C9H7N2S	31257.5
202.0428	C10H8N3S	C10H8N3S	15440.1

Figure 2. MS/MS spectrum of thiabendazole from orange peel with ASAP and an Agilent Q-TOF. The precursor ion peak at 202.0428 m/z and its fragments at 175.0322 m/z and 131.0606 m/z are shown. Mass accuracies are within 3 ppm. The data acquisition time of 0.16 min is circled in green in the top portion of the chromatogram.

Summary

ASAP can be used to rapidly analyze complex samples such as an orange peel directly without sample preparation. This ion source is fully compatible with Agilent TOF and Q-TOF instruments. Typical mass accuracies for MS and MS/MS data are within 3 ppm.

Learn more:

www.agilent.com/chem/qqq

Find a customer center in your country:

www.agilent.com/chem/contactus

U.S. and Canada

1-800-227-9770
agilent_inquiries@agilent.com

Europe

info_agilent@agilent.com

Asia Pacific

inquiry_lsca@agilent.com

For research purposes only and not for use in diagnostic procedures. The information described here is intended for reference and research purposes only. Agilent Technologies offers no guarantee as to the quality and suitability of this data for your specific application.

Information, descriptions and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc. 2010
 Published in the U.S.A. October 25, 2010
 5990-6637EN