

Quick Reference Guide

Installing the Exhaust Fitting (TS600 Module)

Required Material

Loctite® 567 Thread Sealant (included in the exhaust fitting kit).


Installing Procedure



Hand-tightened only.

Do not use tools to tighten any part of the exhaust fitting assembly.

1. Apply a small amount of Loctite 567 to the threaded end of the exhaust fitting hose.
2. Insert the exhaust fitting hose into the module exhaust port.
3. Turn the exhaust fitting clockwise until it is fully seated.
4. Slide the plastic exhaust tube over the exhaust fitting hose.
5. Install the O-ring as shown.
6. Apply a small amount of Loctite 567 at the exposed threaded end of the exhaust fitting hose

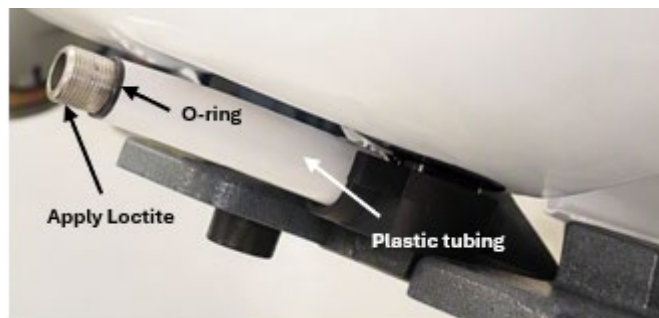

7. and tighten the exhaust flange onto the exhaust fitting hose until secure.

