


Application 125-00


Agilent Fast Natural Gas Analyzer

Technical Overview

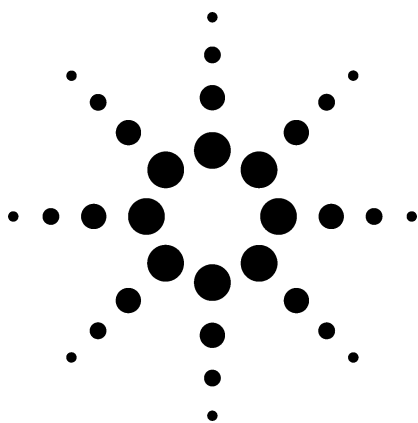

Application Highlights

- A thermal conductivity detector (TCD) to identify air composite (oxygen, nitrogen, and carbon monoxide), methane, carbon dioxide, ethane, propane, isobutane, n-butane, neopentane, isopentane, and n-pentane with an initial C6+ composite backflush to detector.
- 200 ppm lower detection limit for all components except those eluting on the tail of a major preceding constituent.
- System compliant with Gas Processors Association methods 2177 and/or 2261.
- Analysis time is approximately 15 minutes.


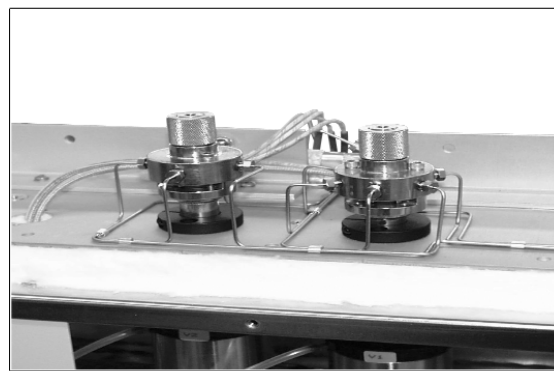


Optional Configurations

- Detailed hydrocarbon analysis of extended natural gas
- TCD/FID/FPD for extended natural gas with trace sulfur analysis
- TCD/FID for extended natural gas with helium or hydrogen

For More Information



For more information on our products and services, visit our Web site at www.agilent.com/chem.

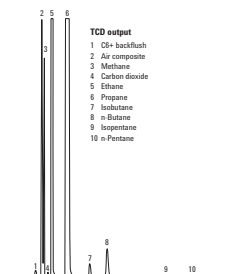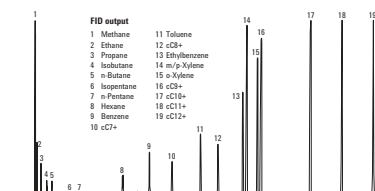

Agilent Technologies



WASSON-ECE
INSTRUMENTATION

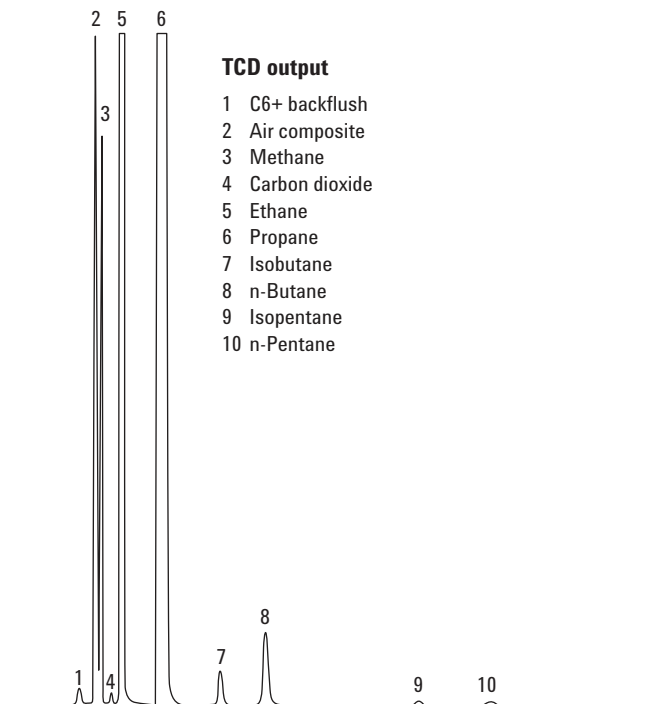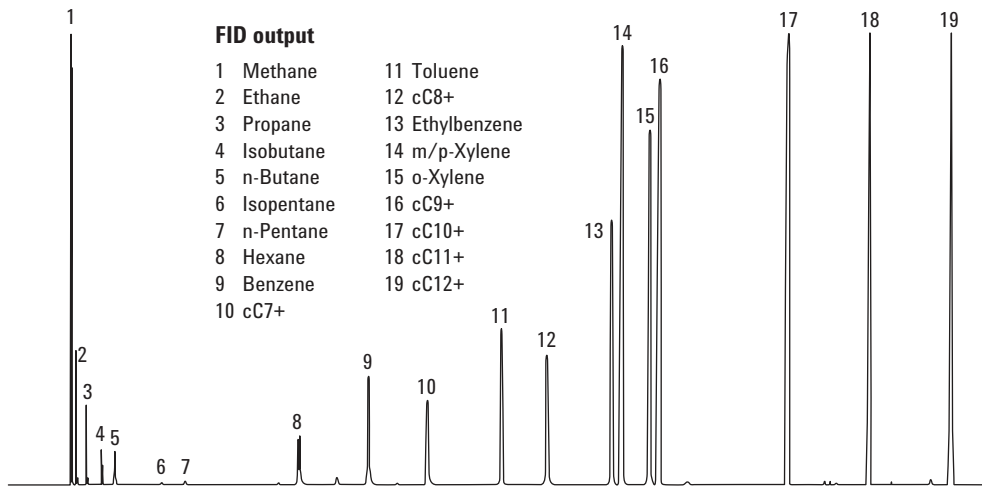

FID and TCD output from Agilent Fast Natural Gas Analyzer.

Agilent shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

Information, descriptions, and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc. 2002

Printed in the USA
December 3, 2002
5988-6725EN

