
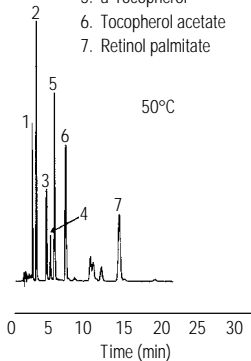


# Fat-Soluble Vitamins on ZORBAX Eclipse XDB-C8

1. Retinol
2. Retinol acetate
3. Vitamin D3
4. g-Tocopherol
5. a-Tocopherol
6. Tocopherol acetate
7. Retinol palmitate

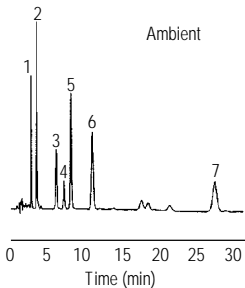

Sample: Fat Soluble Vitamins  
Column: ZORBAX Eclipse XDB-C8, 4.6 x 150 mm, 5  $\mu$ m (P/N: 993967-906)  
Mobile phase: 5/95 Water/MeOH  
Flow rate: 1.0 mL/min  
Temperature: Ambient, 50°C  
Detection: UV 280 nm