


Fusel Oil Simple Standard

DB-624
Custom column
25 m x .53 mm I.D. x 3.0 μ m

Inlet: 250°C, split
FID: 300°C
Carrier: H₂, 50 cm/sec
Oven: 40°C for 5 min.
10°C/min to 250°C

1. acetaldehyde
2. methanol
3. 3-methyl-butanol (isoamyl alcohol)
4. 2-methyl-butanol (active amyl alcohol)

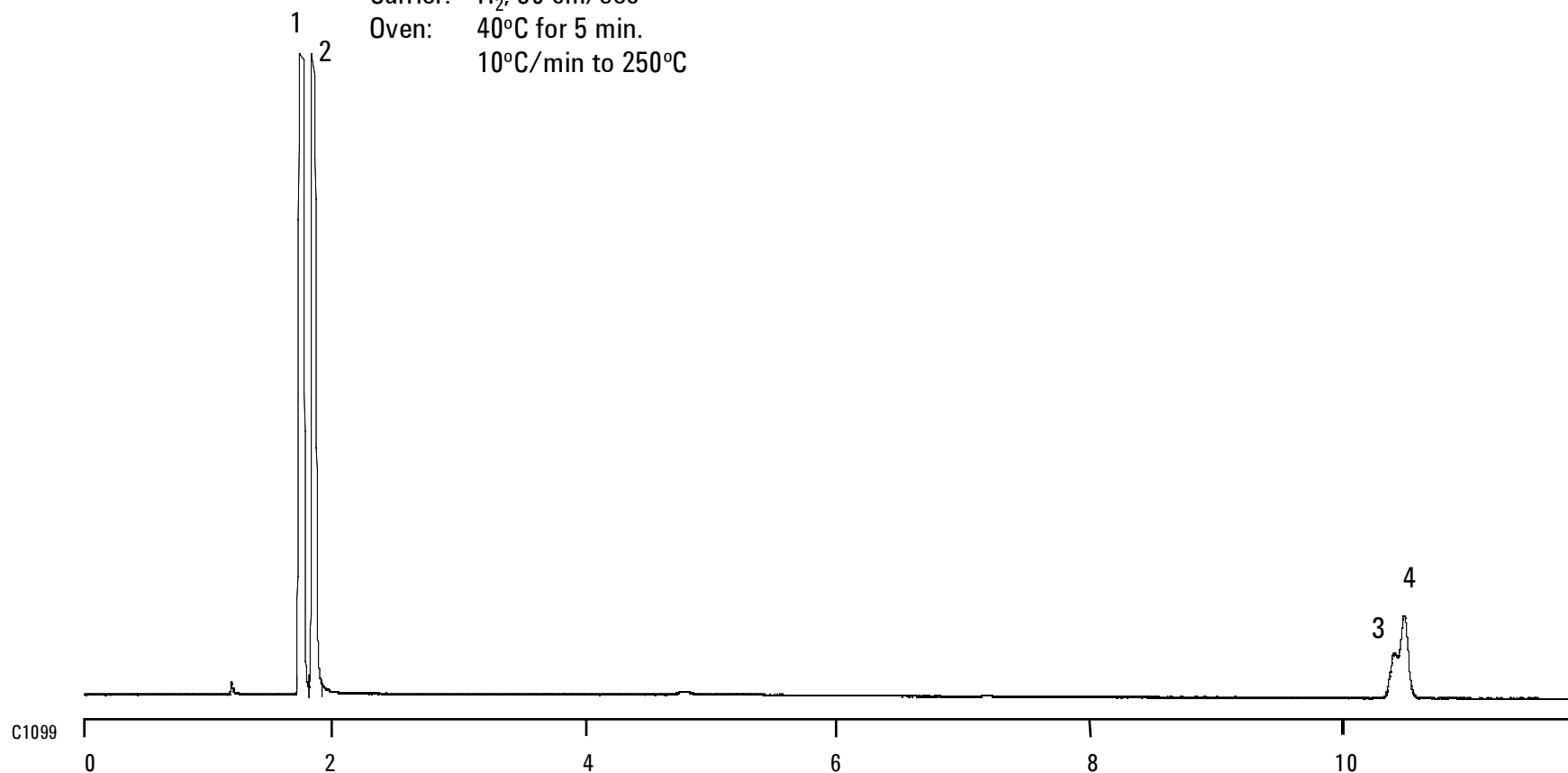

C1099

B-0457

Agilent Restricted



Agilent Technologies

Page 5