


# Fusel Oil Simple Standard

**DB-MTBE**  
**124-0034**  
30 m x .45 mm I.D. x 2.55  $\mu$ m

Inlet: 250°C, split  
FID: 300°C  
Carrier: H<sub>2</sub>, 50 cm/sec  
Oven: 40°C for 5 min.  
10°C/min to 250°C

1. acetaldehyde
2. methanol
3. 3-methyl-butanol (isoamyl alcohol)
4. 2-methyl-butanol (active amyl alcohol)

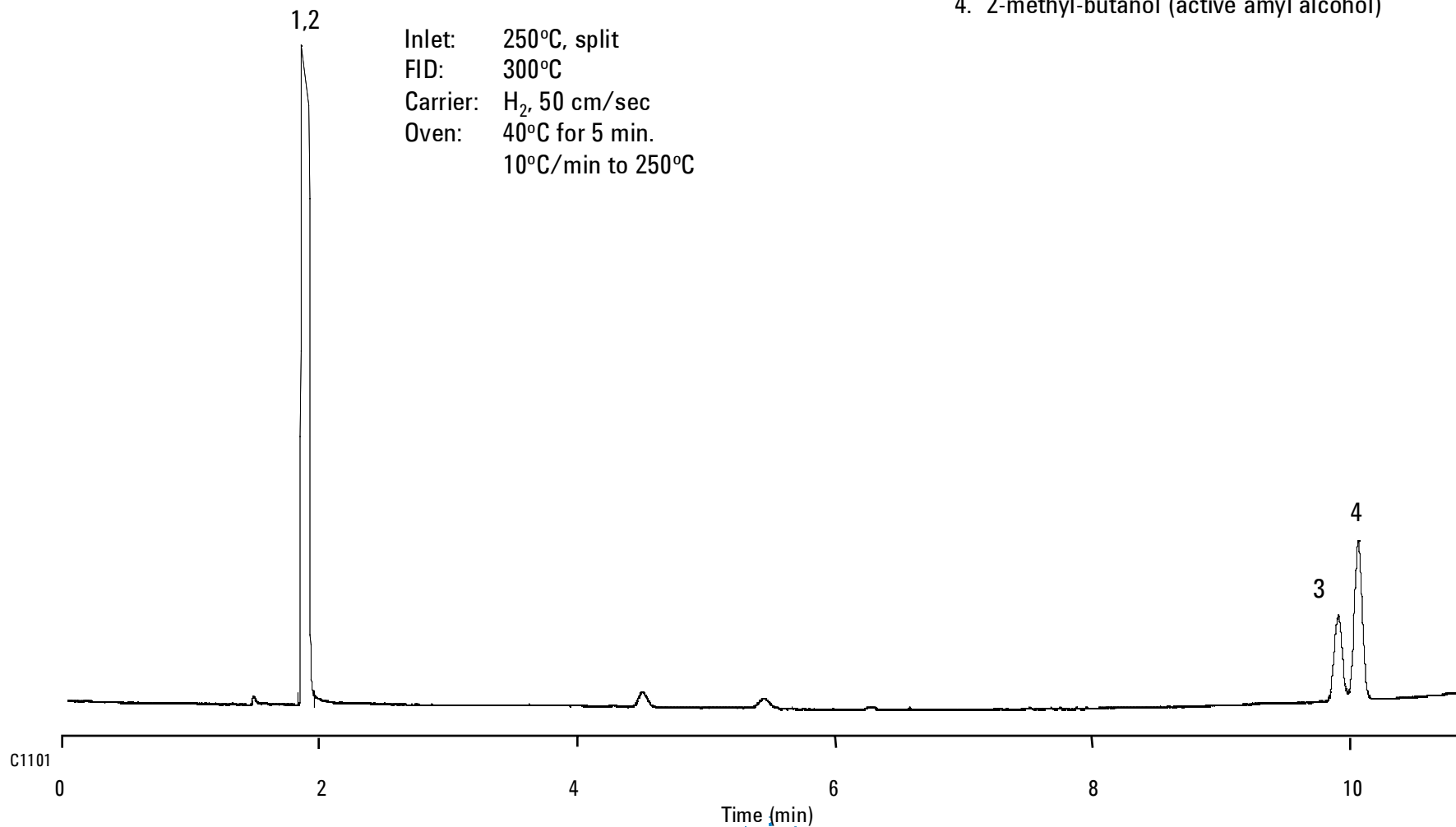

C1101

B-0457

Agilent Restricted



**Agilent Technologies**

Page 2