

## Turpentine Oil

Column: DB-WAX

P/N: 122-7033


Carrier: Hydrogen at 45 cm/  
measured at 70°C

Oven: 70-200°C at 3°/min

Injector: Split 1:150, 220°C

Detector: FID, 220°C

1.  $\alpha$ -Pinene
2. Camphene
3.  $\beta$ -Pinene
4. 3-Carene
5.  $\alpha$ -Phellandrene
6.  $\alpha$ -Terpinene
7. Limonene
8.  $\beta$ -Phellandrene
9.  $\gamma$ -Terpinene
10. *p*-Cymene
11. Terpinolene
12. Caryophyllene
13. Terpinen-4-ol
14.  $\alpha$ -Terpineol

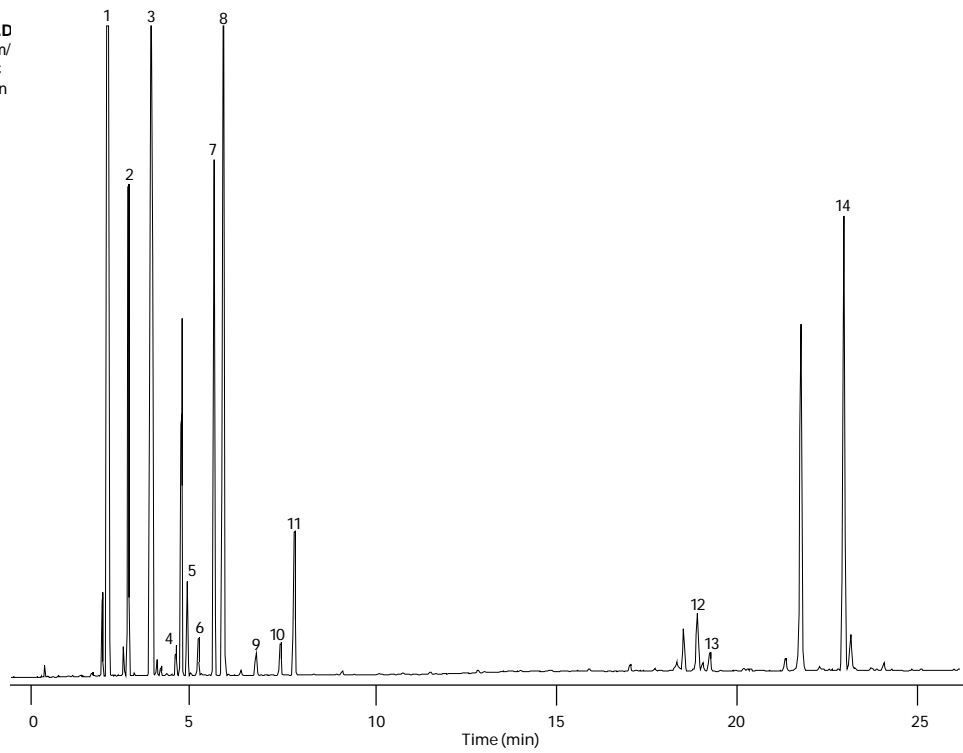

C610