

Turpentine Oil

Column: DB-WAX

P/N: 122-7033


Carrier: Hydrogen at 45 cm/
measured at 70°C

Oven: 70-200°C at 3°/min

Injector: Split 1:150, 220°C

Detector: FID, 220°C

1. α -Pinene
2. Camphene
3. β -Pinene
4. 3-Carene
5. α -Phellandrene
6. α -Terpinene
7. Limonene
8. β -Phellandrene
9. γ -Terpinene
10. *p*-Cymene
11. Terpinolene
12. Caryophyllene
13. Terpinen-4-ol
14. α -Terpineol

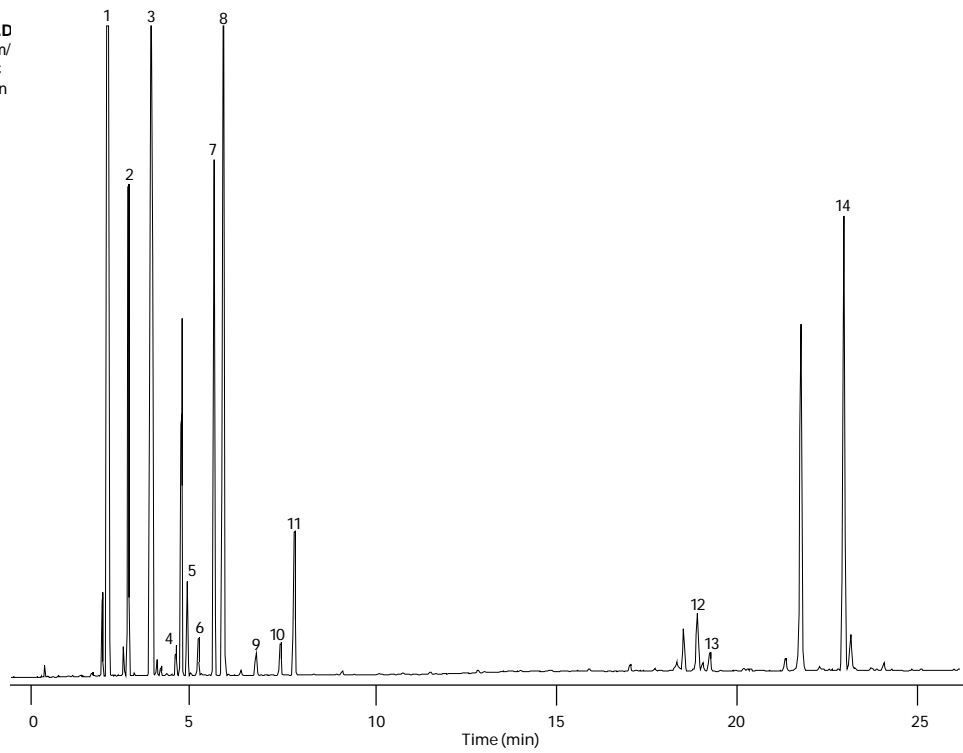

C610