


Alditol Acetates

Column: DB-225
30 m x 0.25 mm I.D., 0.15 μ m
J&W P/N: 122-2231
Carrier: Hydrogen at 36.5 cm/sec
Oven: 220°C isothermal
Injector: Split 1:50, 1 μ L, 225°C
Detector: FID, 250°C
Nitrogen makeup gas at 30 mL/min

1. Rhamnitol
2. Fucitol
3. Ribitol
4. Arabinitol
5. Mannitol
6. Galactitol
7. Glucitol
8. Inositol

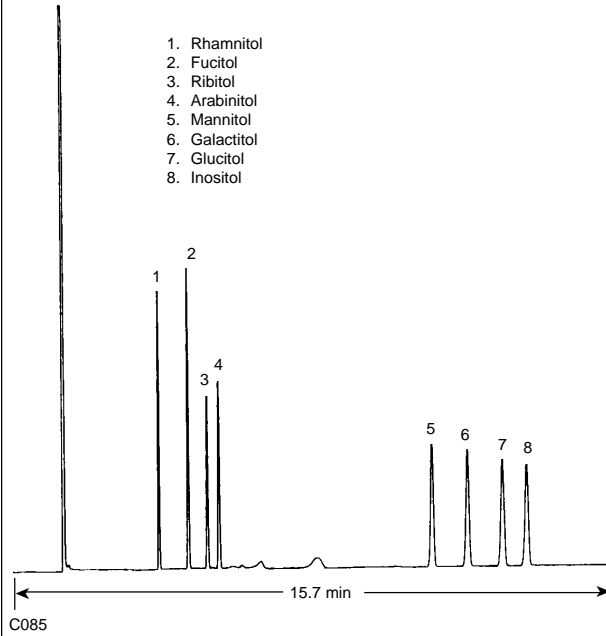

C085