

Anabolic Steroids

Column: DB-1

30 m x 0.25 mm I.D., 0.1 µm

J&W P/N: 122-1031

Carrier: Helium at 40 cm/sec, measured at 180°C


Oven: 180-320°C at 10°/min

320°C for 4 min

Injector: Split 1:40, 2 µL of 0.125 µg/µL each in methanol

Detector: FID, Nitrogen makeup gas at 30 mL/min

1. Dehydroisoandrosterone (Prasterone)
2. 5α-Androstan-17α-ol-3-one (Stanolone)
3. 19-Nortestosterone (Nandrolone)
4. Mesterolone
5. Testosterone
6. 1-Dehydrotestosterone (Boldenone)
7. 17α-Methyltestosterone
8. 1-Dehydro-17-α-methyltestosterone (Methandrostenolone)
9. Norethandrolone
10. 1-Dehydrotestosterone acetate
11. Oxymetholone
12. 19-Nortestosterone-17-propionate
13. 4-Chlortestosterone-17-acetate (Clostebol)
14. Stanozolol
15. 1-Dehydrotestosterone benzoate
16. 19-Nortestosterone-17-decanoate
17. 1-Dehydrotestosterone undecylenate

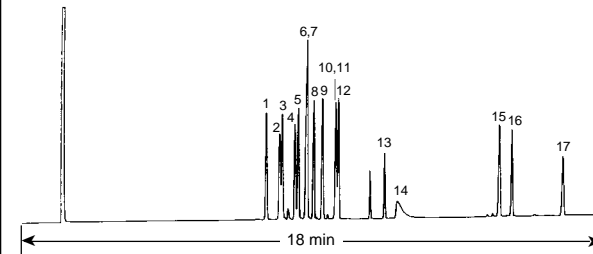

C147