


Triazine Herbicides

Column: DB-17
30 m x 0.25 mm I.D., 0.25 μ m
J&W P/N: 122-1732
Carrier: Helium at 35 cm/sec, measured at 50°C
Oven: 50°C for 1 min
50-100°C at 25°/min
100-250°C at 5°/min
Injector: Splitless, 250°C
40 sec purge activation time
1 μ L of 0.1 μ g/ μ L in methanol
Detector: FID, 300°C
Nitrogen makeup gas at 30 mL/min

1. 2-Nitroxylyene
2. Atraton
3. Propazine
4. Prometon
5. Atrazine
6. Simazine
7. Terbutylazine
8. Secbumeton
9. Prometryne
10. Ametryne
11. Simetryne
12. Terbutryne

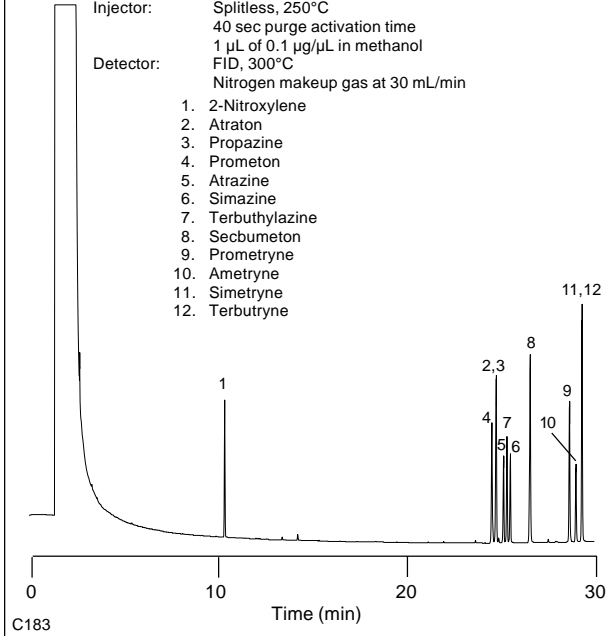

C183