


Triazine Herbicides


Column: DB-35ms
30 m x 0.25 mm I.D., 0.25 µm film
J&W P/N: 122-3832
Carrier: Helium at 33 cm/sec,
measured at 150°C
Oven: 60°C for 1 min
60-340°C at 15°/min
Detector: MSD, 340°C transfer line
Full scan at m/z 45-450
Injector: Hot On-column, 250°C
1 µL of a 10 ng/µL
per component standard

1. Propazine
2. Atrazine
3. Simazine

**Standard 35%
Phenyl**

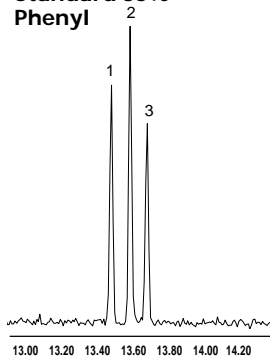

DB-35ms

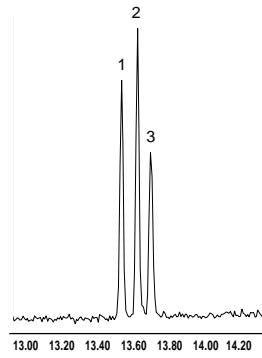

C360