


Narcotics

Column: DB-35ms
30 m x 0.25 mm I.D., 0.25 µm
J&W P/N: 122-3832
Carrier: Helium at 31 cm/sec, measured at 50°C
Oven: 50°C for 0.5 min
50-150°C at 25°/min
150-340°C at 10°/min
Injector: Splitless, 250°C
30 sec purge activation time
Detector: MSD, 300°C transfer line
full scan at m/z 40-380

1. Dextromethorphan
2. Codeine
3. Dihydrocodeine
4. Norcodeine
5. Ethylmorphine
6. Morphine
7. Normorphine
8. 6-Acetylcodeine
9. 6-Monoacetylmorphine
10. Heroin

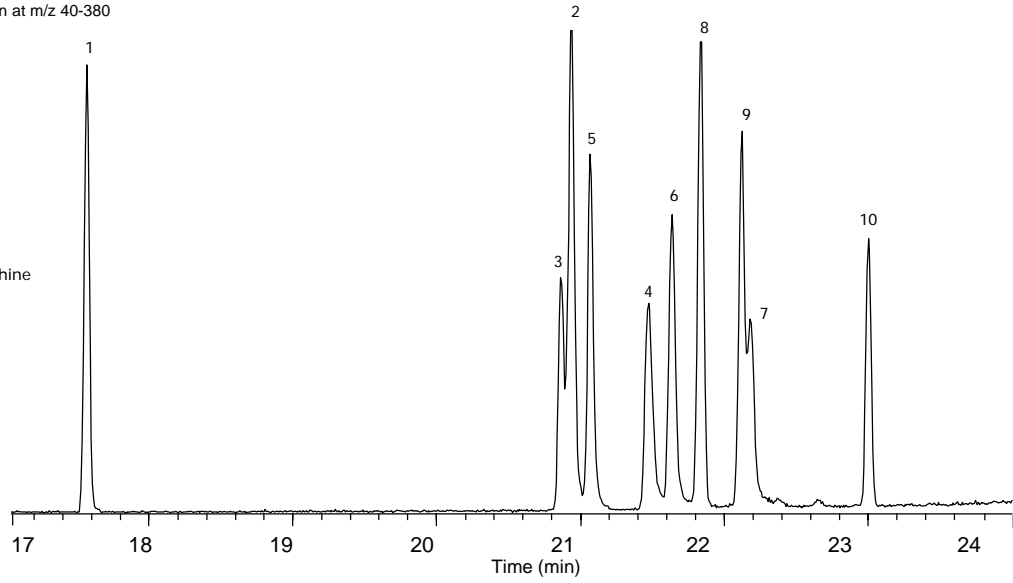

C425