


## Narcotics

**Column:** DB-35ms  
30 m x 0.25 mm I.D., 0.25 µm  
**J&W P/N:** 122-3832  
**Carrier:** Helium at 31 cm/sec, measured at 50°C  
**Oven:** 50°C for 0.5 min  
50-150°C at 25°/min  
150-340°C at 10°/min  
**Injector:** Splitless, 250°C  
30 sec purge activation time  
**Detector:** MSD, 300°C transfer line  
full scan at m/z 40-380

1. Dextromethorphan
2. Codeine
3. Dihydrocodeine
4. Norcodeine
5. Ethylmorphine
6. Morphine
7. Normorphine
8. 6-Acetylcodeine
9. 6-Monoacetylmorphine
10. Heroin

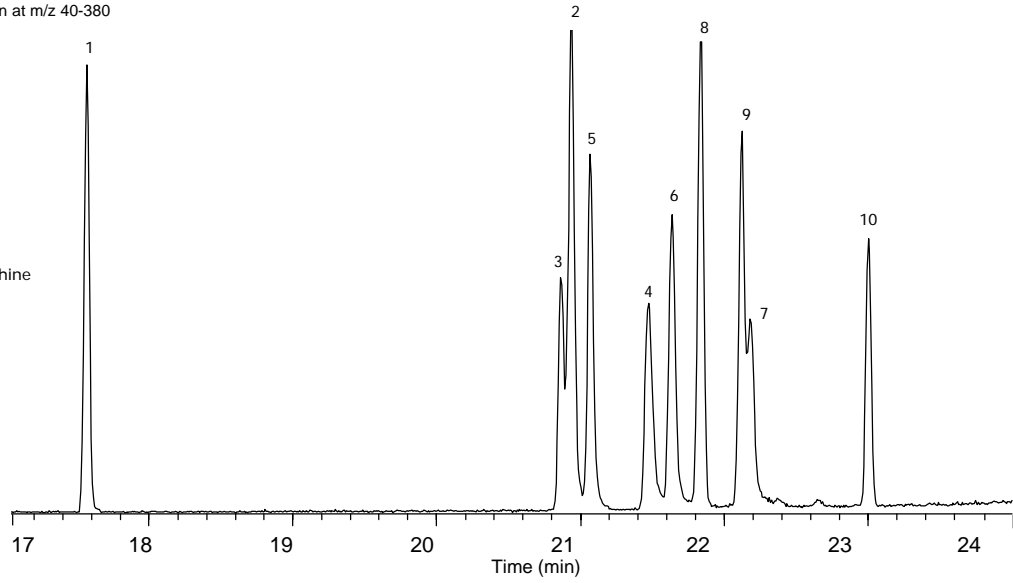

C425