


## Solvents

**Column:** GS-CarbonPLOT  
**30 m x 0.32 mm I.D., 1.5 µm**  
**J&W P/N:** 113-3132  
**Carrier:** Helium at 30 cm/sec (measured at 125°C)  
**Oven:** 125°C for 2 min  
125-350°C at 20°/min  
**Injector:** Split 1:20, 250°C  
**Detector:** MS, SIM mode  
280°C transfer line

1. Methylene chloride
2. Acetone
3. Carbon disulfide
4. Ethyl ether
5. Ethyl acetate
6. Benzene

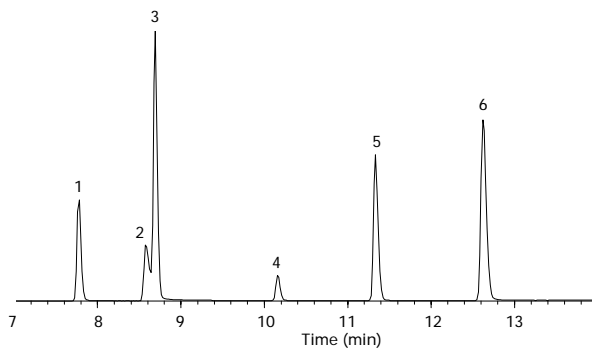

C605