

CLP Pesticides


Column: DB-5
15 m x 0.32 mm I.D., 0.25 µm
10 m DuraGuard Column

Carrier: Helium at 37.2 cm/sec (measured at 50°C)
50°C for 1 min
50-150°C at 25°/min
150-230°C at 4°/min

Injector: Splitless, 250°C
0.5 min purge activation time

Detector: FID, 300°C

- | | |
|-----------------------|------------------------|
| 1. α-BHC | 11. Dieldrin |
| 2. β-BHC | 12. p,p'-DDE |
| 3. γ-BHC | 13. Endrin |
| 4. δ-BHC | 14. Endosulfan II |
| 5. Heptachlor | 15. p,p'-DDT |
| 6. Aldrin | 16. Endrin aldehyde |
| 7. Heptachlor epoxide | 17. Endosulfan sulfate |
| 8. γ-Chlordane | 18. p,p'-DDT |
| 9. Endosulfan I | 19. Endrin ketone |
| 10. α-Chlordane | 20. Methoxychlor |

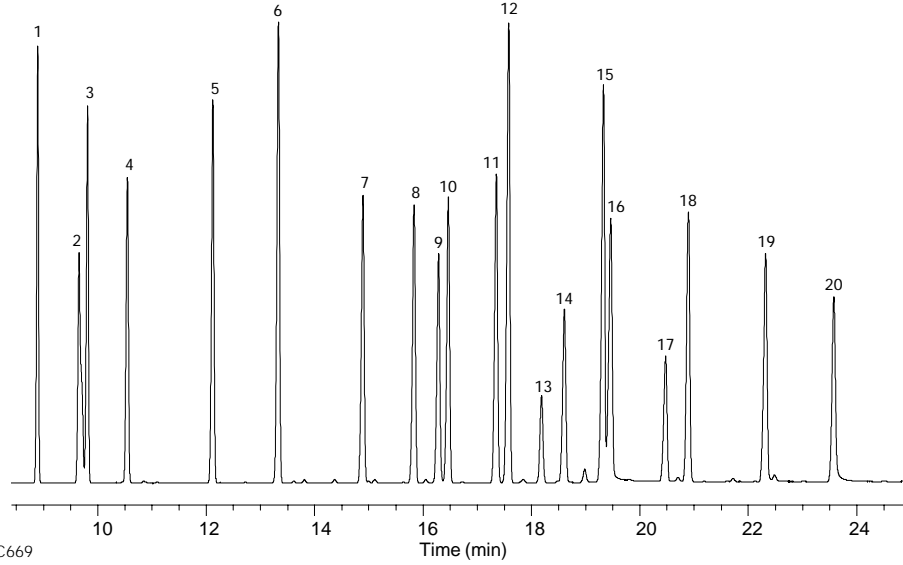

C669

Time (min)