

Acrylate Impurities on DB-PS1701

Column: DB-PS1701

30 m x 0.53 mm, 1.0 µm

J&W P/N: 145-0732

Carrier: Helium, 36.8 cm/sec at 35°C


Oven: 35°C for 5 minutes

35-200°C at 10°/min

Injector: Split, 1:10, 230°C

Detector: FID, 250°C

1. Methanol
2. Ethanol
3. Methyl acetate
4. Ethyl acetate
5. Methyl acrylate
6. Methyl propionate
7. Isobutanol
8. Butanol
9. Ethyl acrylate
10. Ethyl propionate
11. Methyl methacrylate
12. Isopropyl acrylate
13. Isobutyl acetate
14. Butyl acetate
15. Isobutyl acrylate
16. Dibutyl ether
17. Isobutyl propionate
18. Butyl acrylate
19. Butyl propionate

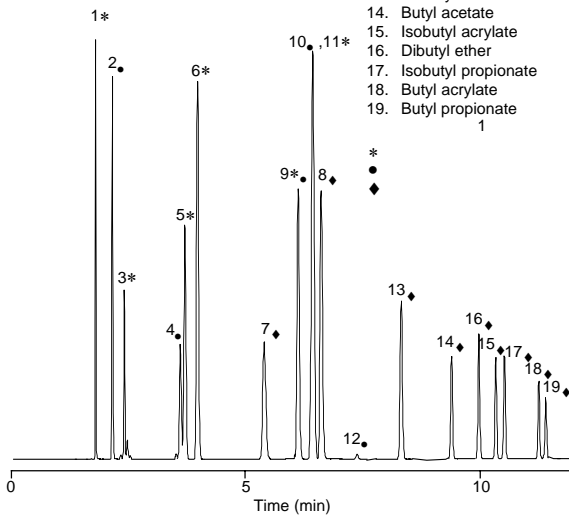

C693