


Extended Sample Compartment

for the Cary 4000, 5000, 6000i, or 7000 UV-Vis-NIR spectrophotometers

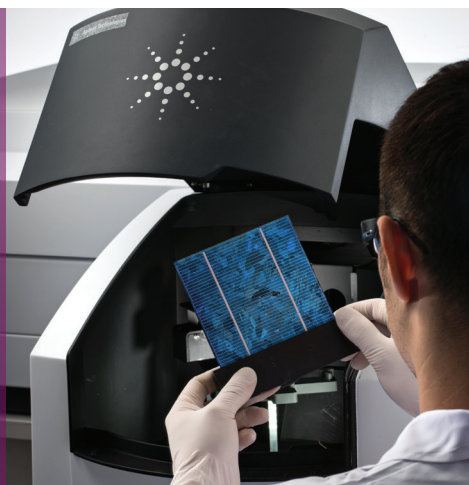


Expanded sampling versatility

An extended sample compartment accessory replaces the standard front wall of the spectrophotometer, increasing the size of the sample compartment. It has excellent light sealing and allows you to measure large samples. The accessory is also required to provide extra sample compartment space to fit some specialized solid/liquid sampling accessories. These include the VW SRA, VASRA, praying mantis, sample transport accessory and internal DRAs.

The extended sample compartment attaches to the front of the spectrophotometer using a slot and grooves. In use, the sample compartment can be accessed by retracting the top sliding lid of the instrument. The instrument sample compartment has internal dimensions of 154 x 150 x 640 mm (W x H x D) when the extended sample compartment is fitted.

Hose fittings are available for use with thermostatted accessories requiring connection to temperature controllers and water baths.

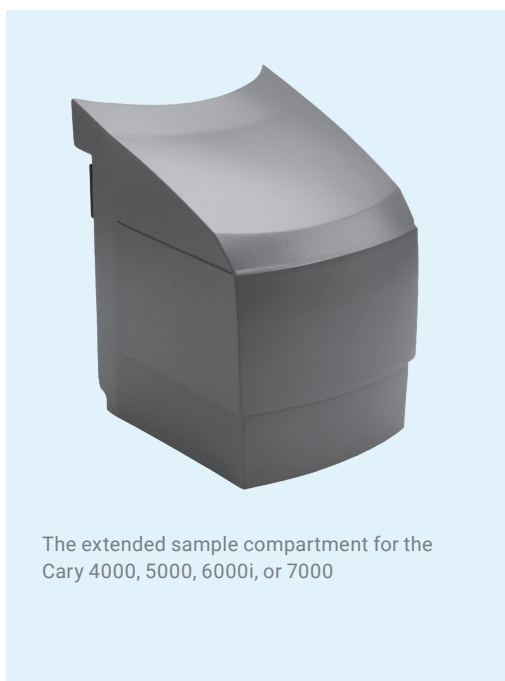

The extended sample compartment for the Cary 4000, 5000, 6000i, or 7000

For more information:

www.agilent.com/chem/uv

DE44462.6517361111

This information is subject to change without notice.

© Agilent Technologies, Inc. 2021
Published in the USA, September 27, 2021
5991-2219EN

